

STERISHEET STERILIZATION WRAPS

STERISHEET 370 SMS



DESCRIPTION

Sterisheet sterilization wraps are best in class Sterile Barrier Systems for CSSDs in hospitals and clinics. Sterisheet SMS products are available as interleaved, bonded, or single sheets.


COMPOSITION

100% synthetic fibers made of polypropylene.

Spunbond / Meltblown / Spunbond

Spunbond network consists of long, strong, and thick polypropylene filaments.

Meltblown layers are made of short and thin polypropylene microfibers.

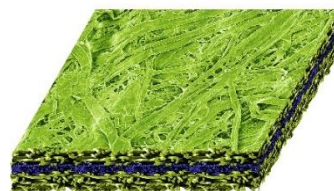

SUITABLE FOR STERILIZATION METHODS

- Steam
- EO
- Vaporized hydrogen peroxide (H₂O₂)

APPLICATION

For large and heavy instrument trays.

SIZES AVAILABLE

Choose standard sizes to optimize your costs.

Other sizes are also available upon request.

75x75 cm 90x90 cm

100x100 cm 120x120 cm

STERISHEET STERILIZATION WRAPS

STERISHEET 370 SMS



COLOURS



PERFORMANCES

Mechanical and bacterial barrier properties are internally tested on a routine basis according to regulatory requirements on all of our products. Combining them with random tests conducted by external and accredited laboratories leads Sterimed to secure the best performances on Sterisheet products.

Mechanical Properties

Preservation of pack integrity from closure till the point of use depends on the material's resistance to tearing, puncturing, and breaching stresses generated all along the distribution with the hospital. Any mechanical weaknesses will increase the risks of event-related ingress of microorganisms into the pack.

Excellent mechanical properties will provide you additional safety while using our materials. Optimal strength and resistance provided in every sheet.

PROPERTIES	TYPICAL
Basis Weight	70 g/m ²
Air permeability	70 l/min/dm ²
Thickness	450 µm
Hydrost Test	50 mBar
Tensile strength MD	3 kN/m
Tensile strength CD	1.3 kN/m
Elongation MD	45 %
Elongation CD	45 %
Burst	325 kPa
Tears MD	9000 mN
Tears CD	14000 mN

Bacterial Properties

Sterilization wraps must prevent microorganisms' ingress inside the package.

To reach this performance, different types of testing have to be performed on the products to reproduce both kinds of ingress vehicles possible:

- Airborne ingress
- Waterborne ingress

For instance, according to the TNO final pack test method, validating the materials & the folding technique, STERISHEET 370 SMS exhibits a BARRIER PROTECTION > 99,99%

STERISHEET STERILIZATION WRAPS

STERISHEET 370 SMS



COMPLIANCE TO STANDARDS

Sterisheet products range is classified as a Class I Medical Device according to the European Medical Device Regulation (MDR 2017/745/EU). Its CE marking illustrates the relevant compliance.

Sterisheet SMS products conform with the standards below:

EN ISO 11607-1

EN 868-2

OUR WRAPS MANUFACTURING CERTIFICATION

ISO 13485 standards

PACKAGING PRIOR TO USE

Sheets presentation is optimized by adjusted folding depending on the size and type of the product. We have carefully tested the best solution to ensure the most convenient handling for end users.

- PRIMARY TRANSPORT PACKAGING

Number of sheets is optimized and wrapped in a transparent polyethylene bag with a quick product ID.

- SECONDARY TRANSPORT PACKAGING

Secondary packaging is a neutral brown color cardboard box providing transportation stress resistance.

LABELING

Product traceability is fully insured through labeling according to regulations on each transport packaging.

STORAGE CONDITIONS

Sterimed recommends the following storage conditions for best performance of sterilization wraps: Storage in a cool, dry location away from direct exposure to natural light, strong artificial light & UV sources.

Cardboard boxes should never be stored in direct contact with the floor. Storage of the products shall be

STERISHEET STERILIZATION WRAPS

STERISHEET 370 SMS



done in areas that are not subject to extreme temperature changes such as in contact with heated objects, vents, or cold walls. As per AAMI ST79 “*Comprehensive guide to steam sterilization and sterility assurance in health care facilities*” recommendations: Before use, hold packaging materials at room temperature (20°C to 23°C) and at a relative humidity ranging from 30% to 60% for a minimum of 2 hours is a good practice for optimum use performances.

USE BY DATE

Provided the above storage conditions are met, the upper limit of the time interval during which the performance characteristics of the sterilization wrap are demonstrated is 5 years of the manufacturing date.

ENVIRONMENTAL IMPACT & WASTE MANAGEMENT

Oil derived product, renewable content 0%.

Disposal as per local regulations after use.